
monzana®

Deuba® Babywaage

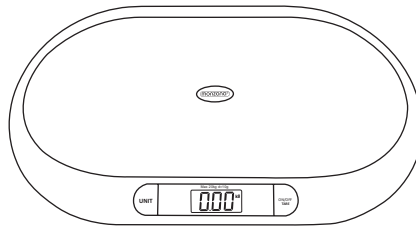Artikel-/Modellnummer:
105021/DBBW001

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INSTRUCTIONS

ATTENTION! Read the safety and assembly instructions carefully, in order to avoid the risk of injury or damaging the product.

Important: Read these instructions carefully and thoroughly. Keep these instructions to read through at a later date. If you give the product to anyone else at some point in the future, please ensure you also pass on this manual

Only for domestic use and not for commercial use!

Keep small parts, plastic bags and other packaging out of the reach of children.

Pay attention especially when unpacking your item. Danger of suffocation!

Check the delivery for completeness. Later complaints cannot be accepted.

Check all elements and parts for damages. Despite careful control, products may be damaged in transit. In this case, do not assemble your product. Defective parts can endanger your health.

Never modify the product! Changes will void the warranty and the product may be unsafe or even dangerous.

Unpack all the parts and elements of the set and place the components on the top of the carton or other clean surface. This protects your new product.

Please also pay attention not to scratch the floor while working with tools.

Make sure that this product is fully assembled before use, as shown in the illustration

For cleaning use a sponge and warm soapy water.

Do not use solvent based-cleaners or detergents - for example bleach. Those can damage the product.

Intended Use


The baby scale is designed exclusively for weighing babies. Use the product only for its designed purpose of use. Manufacturer will not assume any responsibility in case of damages cause of an unappropriated use. Any modification applied by customer on the original product will have a negative impact on its safety, cause danger and warranty will be canceled.

Safety instructions

- Protect your scale from impact with hard objects, moisture, dust, chemicals, toiletries, liquid cosmetics, great temperature fluctuation and closeness to sources of heat (open fires, radiators).
- Do not place any objects on the scale when not in use.
- Strong electromagnetic fields (e.g. cell phones) may impair the

accuracy of the scale.

- Do not overload! The scales can accept a maximum load of 44 lb (20 kg). The weight measurement is displayed in 5 g (0.02 lb) divisions.
- Remove transport safety catch.
- **ATTENTION!** Babies on the scales are to be kept under constant observation in order to avoid falling out..

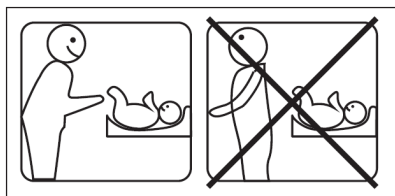Operation

Put the scale on a clean and flat surface. There is the danger that the baby of the scales could fall if you put the scales on a soft or uneven surface. It is important for ensuring that the scale weighs correctly..

Key function

- ON/OFF/TARE key: Power on, Power off, Power tare
- UNIT key:
 - **Select the measuring unit for kg/lb.**
 - **Hold function to fix the weight**

Battery

Use 2pcs AAA battery. Take out the scale, and open up the battery cover, install the batteries to the battery case accordingly. Close the cover.

Weighting step

- 1) Under the state press the key of "ON/OFF/TARE" to power on, the scale will be zero "0.00" appear on the display after 2 seconds.
- 2) Press the UNIT key, you can change "kg", "lb".
- 3) Place the baby on the weighing area. Do this very carefully; do not let the baby fall and take particular care with the baby's head.
- 4) Under the weight state, put the baby on the scale, LCD will show the weight of baby. Do not touch the baby, as this will affect the weight measurement.

- 5) Then push the unit key, to fix the weight., after check the weight, you can push the unit again to “ 0”, continue the next weighing
- 6) Press the key of “ON/OFF/TARE”, LCD shows “0”. Clear the weight of the container, then put the baby into container, the weight will appear on the display.
- 7) Press the key of “ON/OFF/TARE” again, power off.

TIPS:

If you do not want to place the baby directly on the weighing surface, but prefer to use a towel, for example, you have to place the towel on the scales BEFORE switching them on. If the towel (and any other items which you want to deduct again during weighing) weighs over 50 g, you can also use the TARA function.

Error display introduction

- When the scale is over loaded, the “O-LD” will appear on the display.
- When the screen shows “ LO ”, it means that the power isn’t enough, please change the batteries as soon as possible.

Care instructions

If the scale surface is dirty, please clean the scale with soft cloth, do not use any chemical detergent and allow water to get inside the scale.

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This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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Copyright by
Deuba GmbH & Co. KG